

## DATASHEET

The Ambu IV Trainer is for training the insertion of cannulas and catheters, the infusion of fluids and the injection of medication, blood sampling and pulse measurement.

### Replaceable Skin And Veins

Replaceable skin and veins are made of natural latex, which features hole closure after penetration, extending the lifetime of the arm.

### Movable Wrist

The movable wrist and rotating arm ensure an absolutely realistic training product, especially when a catheter is placed on the back of the hand.

### Indicator Panel

The indicator panel which represents each blood vessel, gives immediate feedback, e.g. if an infusion has been correctly accomplished.

### Valve System

The valve system reduces the formation of air pockets in blood vessels.

# Ambu® IV Trainer



## KEY BENEFITS

Replaceable skin and veins

Movable wrist and rotating arm

Indicator panel

Multi purpose carrying case

Fits all AmbuMan® Manikins

Valve system

# Ambu

## SPECIFICATIONS


Product	Dimension
Weight	3 kg
Length	68 cm
Reservoir capacity	260 ml
Carry case dimensions	75 x 15 x 11 cm

## MATERIALS

Product	Material
Skin	Natural latex
Arm and hand	PVC, hard
Veins and artery	Natural latex
Indicator	Silicone
Storage bag	Nylon

## ORDERING SPECIFICATIONS

Item No.	Description
255 001 000	Ambu IV Trainer, complete with infusion stand, red dye concentrate, talcum powder and 5 ml and 50 ml syringes. Supplied in a carry case/washable underlay
<b>Accessories</b>	
255 000 517	Ambu IV Trainer Spare Latex Skin
255 000 501	Ambu IV Trainer Spare Veins, Artery Tubes

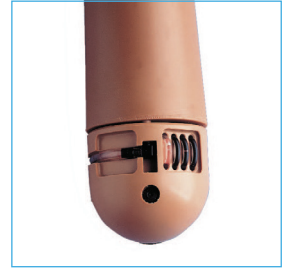

Ambu IV Trainer



Accessories for Ambu IV Trainer



Packed in a robust bag



Ambu, Inc.  
6721 Columbia Gateway Drive  
Suite 200  
Columbia, MD 21046  
Tel. 800 262 8462  
Fax 800 262 8673  
ambuUSA.com

Questions? Contact your Ambu representative  
at **800-262-8462** or visit **ambuUSA.com**



LBL - 008897 - V01 - 01/2023 - Ambu USA. Technical data may be modified without further notice.